



So you can enjoy the important things in life

InnoSpire Mini compressor nebulizer system

InnoSpire Mini allows your patients to nebulize their medications whenever they want and wherever they may be. Compact, lightweight and reliable, it provides independence for people who want to get the most out of life.

An ideal traveling companion, InnoSpire Mini's compact size allows for easy transport while a breath-enhanced reusable nebulizer makes it easy to use. The SideStream Plus reusable nebulizer provides fast, effective treatment and can be used to effectively deliver most commonly prescribed aerosol medications.

PHILIPS
RESPIRONICS

Choose the time and place that's right for you

An optional rechargeable lithium ion battery operates the InnoSpire Mini compressor for up to 90 minutes between charges (and can be fully recharged in approximately 6 hours). A 12V DC (car) adapter enables the compressor to be powered in autos, trucks, boats and recreational vehicles. The compressor can also be powered by the traditional 110V AC wall outlet. InnoSpire Mini is compliant with selected FAA environmental RTCA/DO-160-F standards for airborne equipment as well as the international basic safety and essential performance standard IEC 60601-1 third edition and its amendments.

Accessories include

- SideStream Plus high-efficiency, reusable nebulizer
- Optional rechargeable lithium ion battery
- 12V DC (car) adapter
- AC adapter
- Carrying case
- Spare filters


SideStream Plus breath-enhanced, reusable nebulizer

The SideStream Plus breath-enhanced, reusable nebulizer features an easy-action inspiratory valve that opens on inspiration to boost medication delivery and closes on exhalation, preserving medication and reducing waste.¹

Significantly less pressure is needed to open the inspiratory valve, making it ideal for weak breathers, such as pediatric and elderly patients.



Inspiratory valve

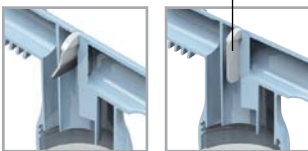

Valve open

1 Performance characterization of the new SideStream Plus breath-enhanced jet nebulizer: J. P. Young, T. J. Hurren, R. H. M. Hatley, T. Dyche. Respironics Respiratory Drug Delivery (UK) Ltd., Chichester, West Sussex, PO20 2FT.

2 Continuous nebulization

3 Free discharge, no nebulizer attached

4 Measured with SideStream Plus nebulizer attached

Product specifications

Nebulization rate ²	0.2 ml. per minute
Nebulization capacity	8 ml.
Mass Median Diameter (MMD) ²	~ 5.0 μ (microns)
Maximum pressure of compressor ³	30 psi
Maximum liter flow of compressor ³	~ 8 liters per minute
Average flow rate of compressor ⁴	3.4 liters per minute @ 12.7 psi operating pressure
Weight	Compressor: 0.93 pounds (0.42 kilograms) Battery: 0.4 pounds (0.188 kilograms)
Size	Compressor: 4.5" x 3.5" x 2.1" Battery: 4.5" x 3.5" x 0.8"
Electrical requirements	12V DC, 1250 mA (compressor) 12V DC, 1700 mAh (battery)
Operating range	41° to 104°F 10 to 95% relative humidity

Ordering information

Product description	Part number	HCPCS code	NDC number
Compressor nebulizer systems			
InnoSpire Mini compressor nebulizer standard	1109289	E0570	08373-9289-00
InnoSpire Mini compressor nebulizer with battery	1109290	E0570	08373-9290-00
Accessories/replacement parts			
InnoSpire Mini rechargeable lithium ion battery	1110726	N/A	08373-0726-00
SideStream Plus high-efficiency reusable nebulizer	HS870	A7005	08373-0870-00
InnoSpire Mini AC power adapter	1109441	N/A	08373-9441-00
InnoSpire Mini 12V DC car adapter	1109455	N/A	08373-9455-00
InnoSpire Mini replacement air filter pack	1109488	A7013	08373-9488-00
Using the InnoSpire Mini compressor with SideStream Plus nebulizer DVD	1110728	N/A	N/A

For custom logo compressors, please contact your Philips Respironics sales representative.

Please visit www.philips.us/innospire-mini for more information

Respironics, InnoSpire Mini and SideStream are trademarks of Respironics, Inc. and its affiliates.

© 2014 Koninklijke Philips N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Geyer JH 4/10/14 MCI 4106138 PN 1112079

Philips Healthcare is part of Royal Philips

Philips Respironics Asia Pacific
+65 6882 5282
Philips Respironics Australia
+61 (2) 9947 0440
Philips Respironics China
+86 400 828 6665
Philips Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
800 345 6443 (toll free, US only)
www.philips.com/respironics

